


B1 小動物室落菌數監測結果

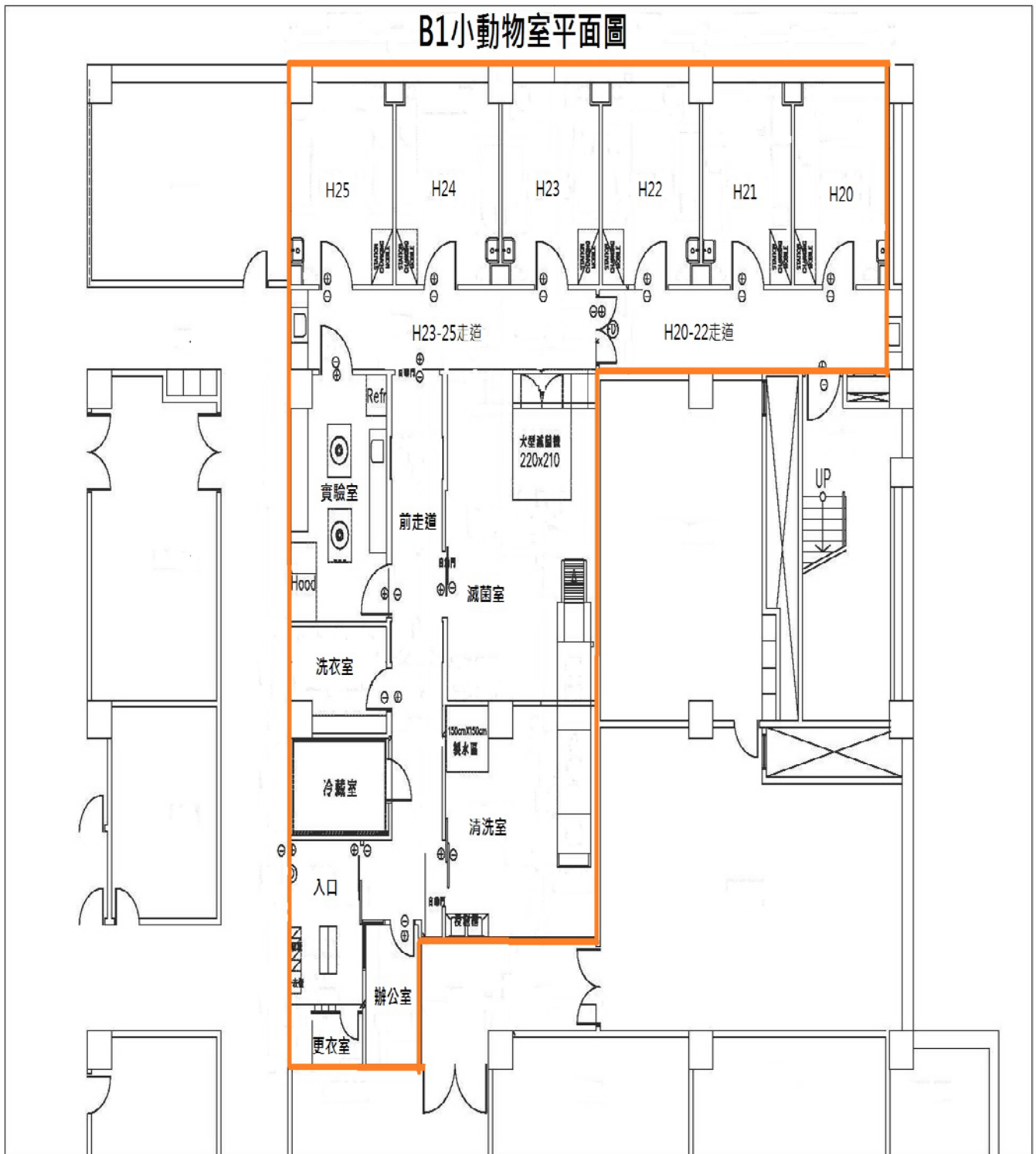

編號	區室	1月	2月	3月	4月	5月	6月	7月	8月	9月	10月	11月	12月
1	H20	0	2	0	0	0	0	0	0	0	0	0	0
2	H21	1	1	0	0	0	0	0	0	0	0	0	2
3	H22	0	0	0	1	0	0	0	0	2	0	0	0
4	H23	1	1	0	1	1	0	1	1	0	0	0	0
5	H24	0	4	1	1	0	2	0	0	0	1	0	0
6	H25	0	2	1	2	0	0	0	1	1	0	0	0
7	20-22 室走道	0	0	2	0	0	1	0	3	2	1	0	0
8	23-25 室走道	1	3	1	1	0	0	1	1	3	0	1	1
9	實驗室	1	1	0	0	0	0	0	2	1	0	1	1
10	洗衣室	0	0	0	0	0	2	1	2	10	2	2	1
11	滅菌室	1	0	0	1	0	0	0	1	2	1	0	1
12	清洗室	0	0	1	1	0	0	2	0	5	2	0	2
13	前走道	0	0	0	1	1	0	2	2	3	0	2	3
14	辦公室	0	2	6	0	1	1	1	12	9	4	0	3
15	更衣區	2	5	3	4	17	0	7	12	7	8	11	10
16	入口	5	4	6	4	11	0	5	15	8	6	5	4
17	20 操作台	1	0	0	0	0	0	0	0	0	0	0	0
18	21 操作台	0	0	1	0	0	0	0	0	0	0	0	0
19	22 操作台	0	0	0	0	0	0	0	0	0	0	0	0
20	23 操作台	0	0	0	0	0	0	0	0	0	0	0	0
21	24 操作台	1	0	0	0	0	0	0	0	0	0	0	0
22	25 操作台	0	0	0	0	0	0	0	0	0	0	0	0

1. BBL™ Stacker™ Plate TSB II , 5% Sheep blood (Becton, Dickinson)。
2. 方法：打開 Plate 上蓋，置於房間中央 30 分鐘，30 分鐘後蓋好上蓋，置入 37°C 培養箱中培養 48 小時，取出後計算落菌數量